



# 4' & 5' Protector® Evidence Drying Cabinets

## SPECIFICATIONS

### ALL MODELS FEATURE:

- Type 304 brushed stainless steel interior with 2" OD floor drain, stainless steel hanging rod and three stainless steel collapsible/removable wire mesh shelves, 8.0" w x 19.0" d (20.3 x 48.3 cm)
- Epoxy-coated steel exterior
- Safety glass double door with magnetic seal that accommodates a user-supplied padlock or tamper-evident tag
- Inlet and exhaust roughing prefilters with MERV 11 rating
- Exhaust HEPA filter, 99.99% efficient on particles 0.3 micron
- Activated organic vapor carbon exhaust filter, 7.5 lbs. (3.4 kg) porous coconut shell media
- Built-in blower with solid state speed control
- Red, filter system alarm light to alert operator to low airflow conditions
- ETL listed (115V and 230V North America models)
- CE conformity marking (230V International models)
- 8' (2.4 meter) power cord with plug
- Shipping weight 585 lbs. (265 kg) for 4' and 675 lbs. (306 kg) for 5'

### MODELS WITH UV LIGHT INCLUDE:

- 254 nm, key-activated UV lamp,
- Storage area, 44.25" w x 26.5" d x 12.0" h (112.4 x 67.3 x 30.5 cm) for 4' models and 56.25" w x 26.5" d x 12.0" h (142.9 x 67.3 x 30.5 cm) for 5' models

### MODELS WITH WASHDOWN INCLUDE:

- Spray nozzle with adjustable stream and 6' (1.8 meter) hose
- 10' (3.5 meter) rubber drain hose preplumbed to drain
- Storage area, 27.5" w x 26.5" d x 12.0" h (69.9 x 67.3 x 30.5 cm) for 4' models and 39.5" w x 26.5" d x 12.0" h (100.3 x 67.3 x 30.5 cm) for 5' models

## ORDERING INFORMATION

Catalog #	4' Model	5' Model	UV Light	Washdown	Electrical Requirements	Plug Type
3404000	X		X		115 volts, 60 Hz, 6 amps	115V, 13 amps
3404001	X		X		230 volts, 50/60 Hz, 3 amps	North America, 230 Volts
3404002*	X		X		230 volts, 50/60 Hz, 3 amps	British (UK)
3404003*	X		X		230 volts, 50/60 Hz, 3 amps	Schuko
3404004*	X		X		230 volts, 50/60 Hz, 3 amps	China/Australia
3404010	X			X	115 volts, 60 Hz, 6 amps	115V, 13 amps
3404011	X			X	230 volts, 50/60 Hz, 3 amps	North America, 230 Volts
3404012*	X			X	230 volts, 50/60 Hz, 3 amps	British (UK)
3404013*	X			X	230 volts, 50/60 Hz, 3 amps	Schuko
3404014*	X			X	230 volts, 50/60 Hz, 3 amps	China/Australia
3405000		X	X		115 volts, 60 Hz, 6 amps	115V, 13 amps
3405001		X	X		230 volts, 50/60 Hz, 3 amps	North America, 230 Volts
3405002*		X	X		230 volts, 50/60 Hz, 3 amps	British (UK)
3405003*		X	X		230 volts, 50/60 Hz, 3 amps	Schuko
3405004*		X	X		230 volts, 50/60 Hz, 3 amps	China/Australia
3405010		X		X	115 volts, 60 Hz, 6 amps	115V, 13 amps
3405011		X		X	230 volts, 50/60 Hz, 3 amps	North America, 230 Volts
3405012*		X		X	230 volts, 50/60 Hz, 3 amps	British (UK)
3405013*		X		X	230 volts, 50/60 Hz, 3 amps	Schuko
3405014*		X		X	230 volts, 50/60 Hz, 3 amps	China/Australia



Shown above: 5' Protector Evidence Drying Cabinet 3405010

Note: See reverse side for ordering information on accessories including paper liners and cleaners.



## 4' & 5' Protector® Evidence Drying Cabinets

### Accessories & Dimensional Data

#### ACCESSORIES

**Starter Kit.** Includes six Inlet Prefilters, six Exhaust Prefilters, 10 Hanging Clips, 10 Security Tags, 10 Paper Liners, one 32 ounce (946 ml) bottle of Germicidal Cleaner.

Catalog Number	Description	Shipping Weight
3869910	For 4' Cabinets	16 lbs. (7.3 kg)
3869920	For 5' Cabinets	17 lbs. (7.7 kg)

**3869600 Exhaust Connection Kit.** Allows the cabinet to be connected to house exhaust or other exhaust system. Includes dry powder epoxy-coated steel connector with 5" (127 mm) diameter collar, 14' (36 cm) of 5" diameter (127 mm) flexible duct, and two clamps. The user-supplied exhaust system must be capable of exhausting 150 CFM (4.2 m³/min) at a vacuum of 0.5" of water (125 Pa). Shipping weight 10 lbs. (5 kg)

**3928200 Drain Kit.** Includes 10' (305 cm) of drain hose and fitting to connect to floor drain. Shipping weight 10 lbs. (5 kg)

**Paper Liners.** White paper is precut to fit the floor of the cabinet to collect material that may fall from evidence. Paper is polyethylene coated on one side.

Catalog Number	Description	Shipping Weight
3869802	For 4' Cabinet. Package of 25.	12 lbs. (5.5 kg)
3894201	For 4' Cabinet. Package of 100.	25 lbs. (11 kg)
3869804	For 5' Cabinet. Package of 25.	12 lbs. (5.5 kg)
3894202	For 5' Cabinet. Package of 100.	25 lbs. (11 kg)



**3902400 Security Tags (left).** Tamper-evident, label tags irreversibly indicate if the cabinet door has been opened during the drying cycle. Package of 100. Shipping weight 1 lb. (0.5 kg)



**3869500 Hanging Clips (right).** Plastic clips allow small articles to be hung from the cabinet's hanger rod or shelf. Package of 10. Shipping weight 1 lb. (0.5 kg)

**Germicidal Cleaners.** Clean, disinfect and deodorize. May be used on hard, non-porous surfaces including stainless steel.

Catalog Number	Description	Towel Size	Shipping Weight
3872400	100 individually packaged single use wipes	6.00" x 8.00"	2 lbs. (1 kg)
3872500	90 wipe dispensing tub	6.00" x 10.50"	2 lbs. (1 kg)
3901800	Bottle, 32 ounces. Includes screw-on spray head dispenser.		2 lbs. (1 kg)
3901700	Container, 1 gallon. For refilling spray bottle.		10 lbs. (5 kg)

#### FILTERS

**6938100 Replacement HEPA Exhaust Filter.** 99.99% efficient on particulates 0.3 micron. One HEPA Exhaust Filter is included with each Protector Evidence Drying Cabinet. Shipping weight 9.0 lbs. (4.0 kg)

**6938200 Replacement Activated Carbon Filter.** Contains 7.5 lbs. (3.4 kg) of porous, coconut shell carbon. One Activated Carbon Filter is included with each Protector Evidence Drying Cabinet. Shipping weight 12.5 lbs. (4.7 kg)

**Replacement Door Prefilters, Inlet.** Polyester prefilters trap larger particulates in air entering the cabinet. MERV 11 rating. Precut to fit the door of the cabinet. Package of 6.

Catalog Number	Description	Shipping Weight
3872601	For 4' Cabinets	1 lb. (0.5 kg)
3872602	For 5' Cabinets	1 lb. (0.5 kg)

**3872700 Replacement Cabinet Prefilters, Exhaust.** Polyester prefilters trap larger particulates in air exhausted from the cabinet to extend the life of the exhaust HEPA filter. MERV 11 rating. Precut to fit the inside top of the cabinet interior. Package of 6. Shipping weight 1 lb. (0.5 kg)



**Single-Use Door Evidence Prefilter.** May be examined and/or retained for trace evidence. Single-use. Polyester. Pre-cut to fit the door prefilter frame. Package of 10. Shipping weight 1 lb. (0.5 kg)

Catalog Number	Description	Shipping Weight
3951822	For 4' Cabinets	1 lb. (0.5 kg)
3951832	For 5' Cabinets	1 lb. (0.5 kg)

**3951803 Single-Use Cabinet Evidence Prefilter.** May be examined and/or retained for trace evidence. Single-use. Polyester. Pre-cut to fit the cabinet prefilter frame. Package of 10. Shipping weight 1 lb. (0.5 kg)

#### DIMENSIONAL DATA

SIZE	"A"	"B"
4'	48.00" (121.9 cm)	44.25" (112.4 cm)
5'	60.00" (152.4 cm)	56.25" (142.9 cm)

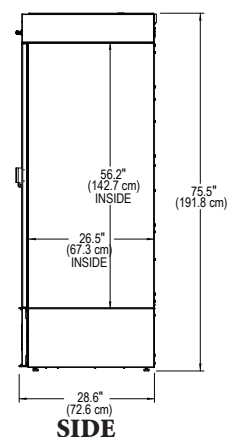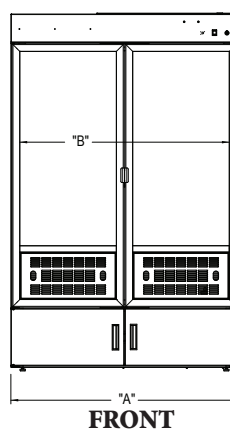

Protecting your  
laboratory environment  
**LABCONCO**

**Labconco Corporation**

8811 Prospect Avenue • Kansas City, Missouri 64132  
816-333-8811 or 800-821-5525 • www.labconco.com

