

GENERAL DESCRIPTION

4' RXPERT DOUBLE FILTERED BALANCE SYSTEM - 115 V

Model: 3971492

Standard Features

- * Deep 30" interior to accommodate large micro and analytical balances
- * Ergonomic air foil with aerodynamic Clean-Sweep® airflow openings
- * Upper dilution air supply
- * Upper containment sash foil
- * Zone-perforated, pivoting rear baffle®
- * Side-entry air foils
- Two exhaust HEPA filter, 99.99% efficient on particles 0.3 micron in size
- * Each HEPA filter 100% independently leak-detection scannable without removing impellers or filters
- True bag-in/bag-out filter disposal system
- Built-in blower with dynamically balanced, vibration-isolated motorized impeller and speed control
- Pressure gauge to monitor filter loading
- Glacier white and gray, dry powder epoxy-coated aluminum frame and steel rear plenum and baffle.
- Pre-wired fluorescent lamp and switches for blower and light.
- Guardian LED Analog Airflow Monitor
- 1/4" thick tempered safety glass front sash and sides
- Two utility ports for pass-through of electrical cords
- Solid Black Epoxy Dished Work Surface


* Labconco Exclusive Feature

CONFORMITY


- * USP < 800 > non-sterile compounding
- * SEFA 1
- * Modified ASHRAE 110
- * ANSI Z9.5
- * UL 61010-1
- * CAN/CSA C22.2 No. 1010.1

SOLD SEPARATELY


- * Supporting Base

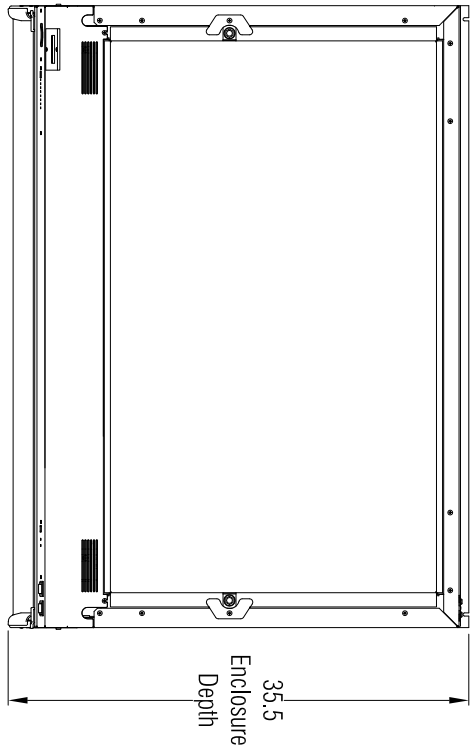

Top View

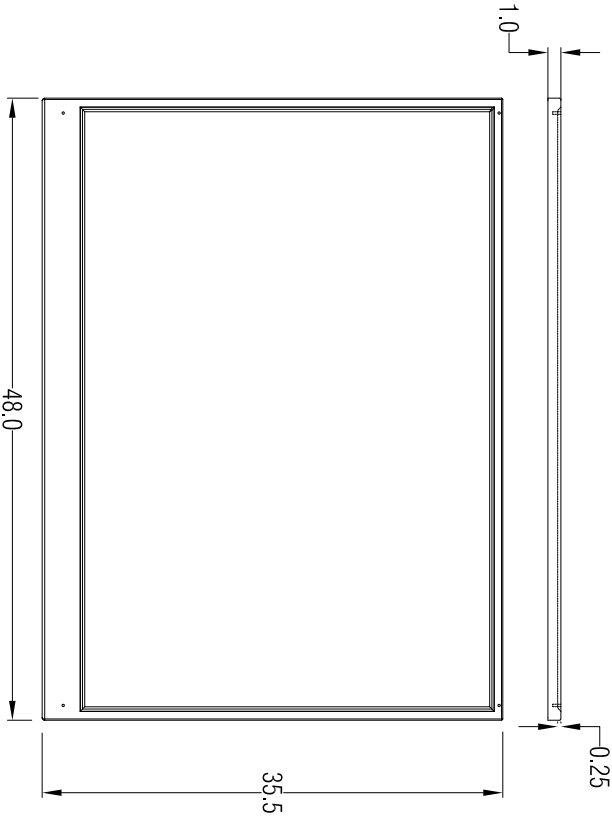

Work Surface

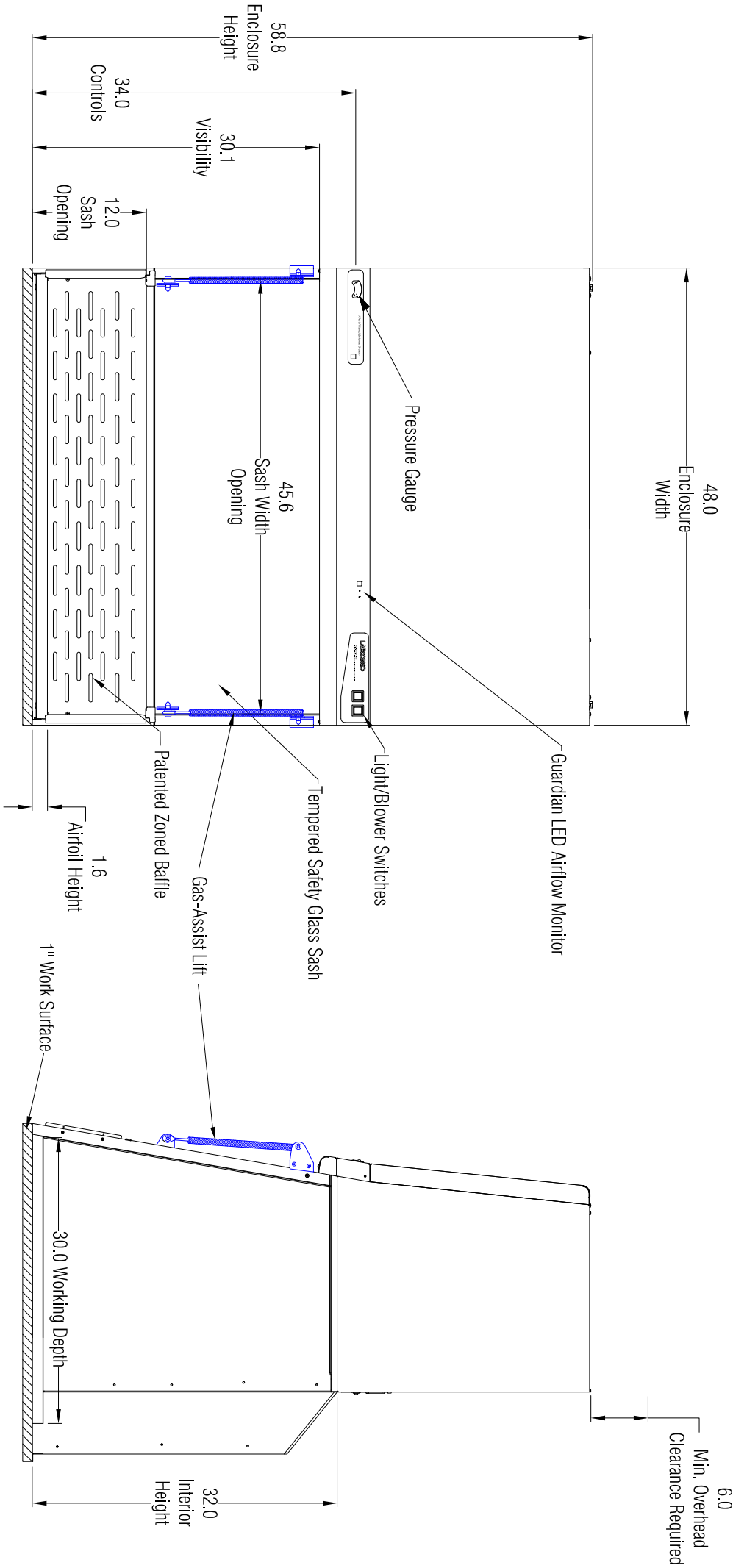

Front View

Side View

JOB NAME:



LABCONCO CORPORATION • KANSAS CITY MO USA 64132
800-821-5525 / 816-333-8811 • labconco.com

REVISIONS:

No. Date

DATE:
12-13-2022

DRAWING:

3971492

DRAWN BY:

C. Alumbough

Order Number:

PROJECT No.

SHEET No.