

GENERAL DESCRIPTION

5' RXPERT DOUBLE FILTERED BALANCE SYSTEM - 115 V
Model: 3971592

Standard Features


- * 30" interior to accommodate large micro and analytical balances
- * Ergonomic air foil with aerodynamic Clean-Sweep® airflow openings
- * Upper dilution air supply
- * Upper containment sash foil
- * Zone-perforated, pivoting rear baffle®
- * Side-entry air foils
- Two exhaust HEPA filter, 99.99% efficient on particles 0.3 micron in size
- * Each HEPA filter 100% independently leak-detection scannable without removing impellers or filters
- True bag-in/bag-out filter disposal system
- Built-in blower with dynmically balanced, vibration-isolated motorized impeller and speed control
- Pressure gauge to monitor filter loading
- Glacier white and gray. dry powder epoxy-coated aluminum frame and steel rear plenum and baffle.
- Pre-wired fluorescent lamp and switches for blower and light.
- Guardian LED Analog Airflow Monitor
- 1/4" thick tempered safety glass front sash and sides
- Two utility ports for pass-through of electrical cords
- Solid Epoxy Work Surface

* Labconco Exclusive Feature

Required Accessories: (Ordered Separately)


* Supporting Base

- CONFORMITY**
- * USP <800> non-sterile compounding
 - * SEFA 1
 - * Modified ASHRAE 110
 - * ANSI Z9.5
 - * UL 61010-1
 - * CAN/CSA C22.2 No. 1010.1

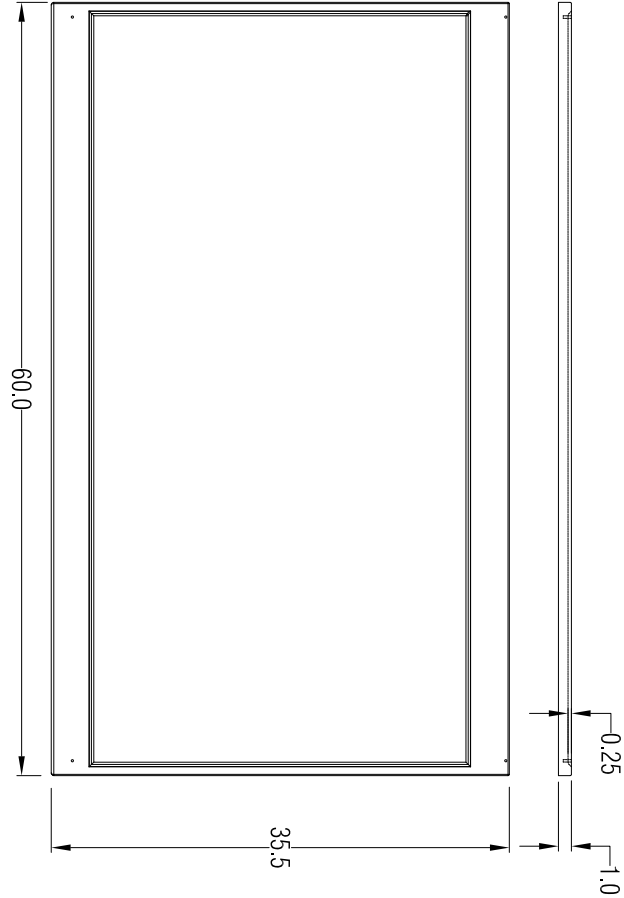


Top View

Solid Epoxy Work Surface

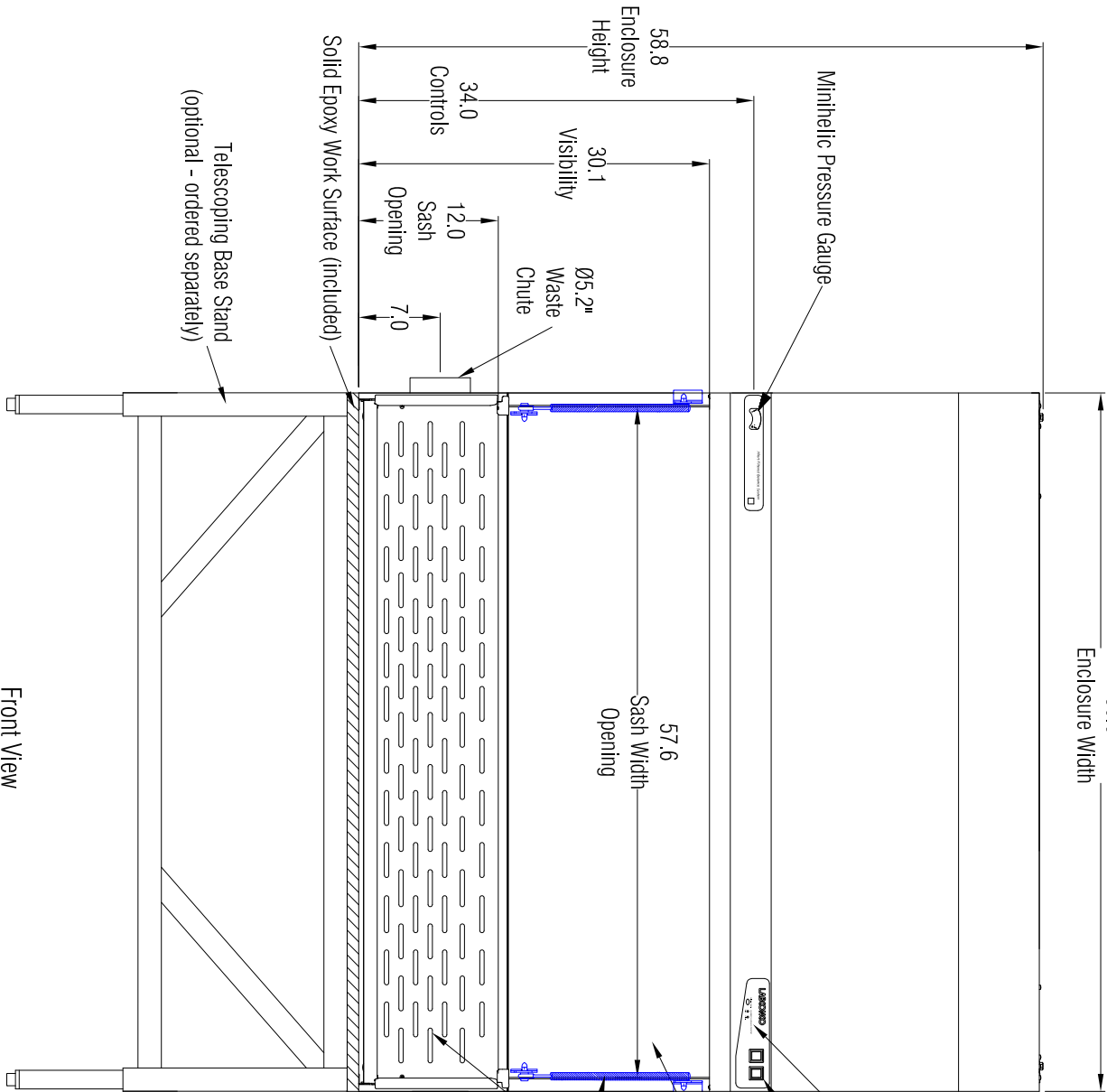

Front View

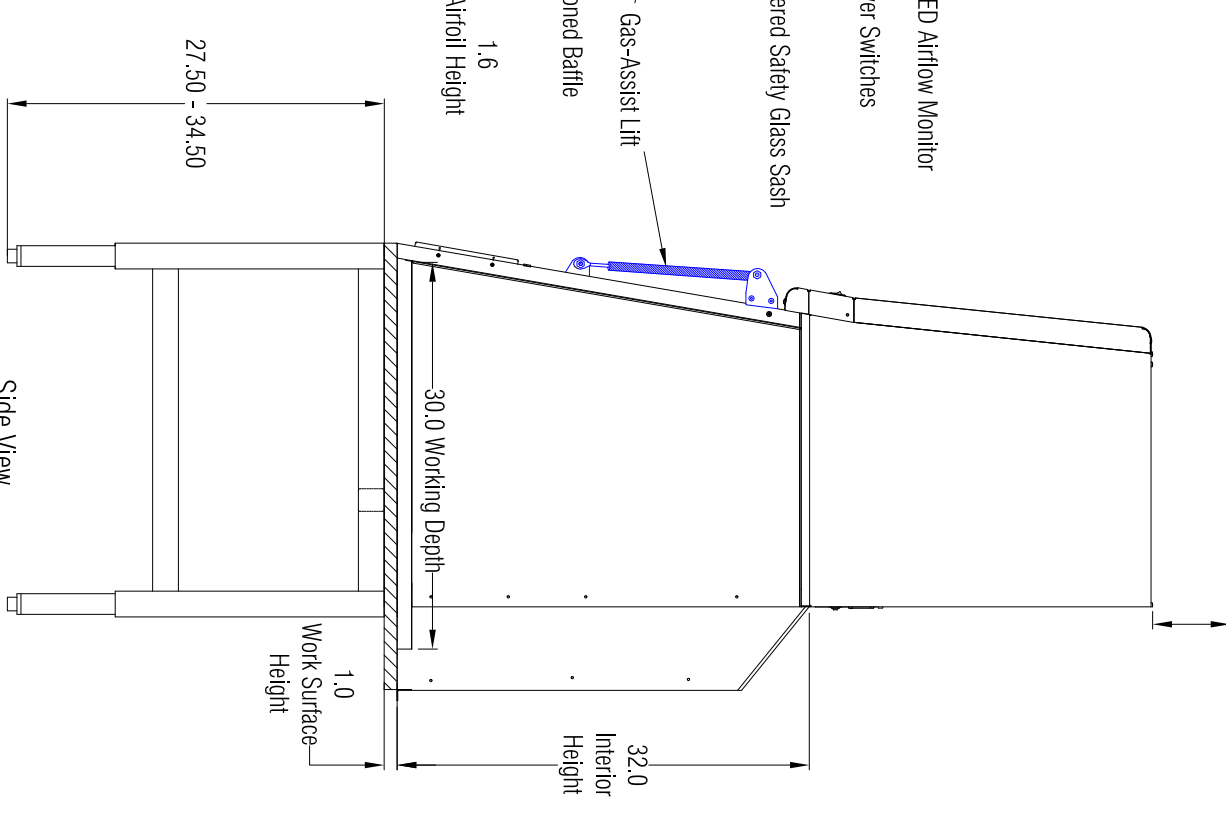

Side View



LABCONCO CORPORATION • KANSAS CITY MO USA 64132
800-821-5525 / 816-333-8811 • labconco.com

JOB NAME:

REVISIONS:	
No.	Date

DATE:
09-29-2023

DRAWING:
3971592

DRAWN BY:
C. Alumbough

Order Number:

PROJECT No.

SHEET No.