



ISO View

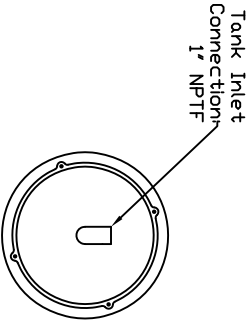

Front View

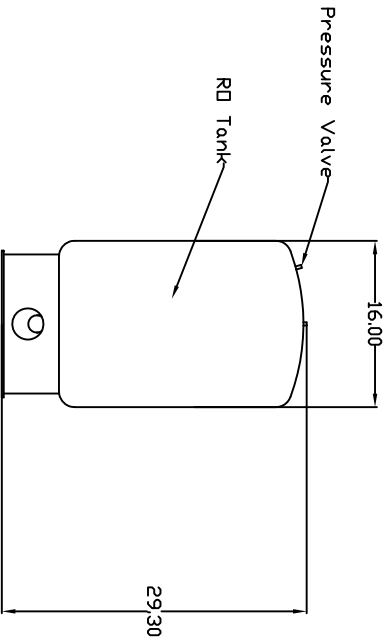

Side View

WaterPro RO Bladder Tank, 20 Gallon  
Model # 9201000

- Provides additional storage of reverse osmosis water produced by the WaterPro RO Station.
- Pressurized water can be distributed to a polishing station or to a glassware washer through either of two pressurized ports.
- Demand of the WaterPro RO Bladder Tank reduces the system pressure, which turns on the WaterPro RO Station to refill the tank.
- If demand is such that the pressurized WaterPro RO Bladder Tank empties, this system can also draw water from the 17L gravity feed storage tank located on the WaterPro RO Station.
- Actual capacity amended due to pressures: 12.87 gal (48.7 L).
- Overall dimensions: 16" dia x 29.3" h
- Shipping weight: 38 lbs (17 kg)
- *Tubing and fittings not shown on drawing.*

Bladder Tank for:



Date:  
01-06-2023

Drawing #  
9201000

Drawn By:  
C. Alumbaugh

SCALE:      Scaled to fit